

How much do target standards improve this calibration transfer?


CalibShift — transfer to changed acquisition conditions | Version 0.5.2

Execution completed. Exploratory results; selected uncertainty measures unavailable. Units: % w/w.


Group RMSE	MAE	Bias	Usable pairs
1.43	1.213	-0.7907	9

Transfer direction is X1 to X2. The selected Bias correction uses 2 target standards and gives group RMSE 1.43 and bias -0.7907. The budget plot uses the common held-out set. Unavailable fits and standard identities are recorded in extended tables. Increasing a budget or fitting a slope can worsen performance. This protocol does not establish transfer across arbitrary instruments or laboratories. Primary score and scatter: 9 paired rows (common comparison cohort); 9 pairs are available for this method in the full export. Group RMSE weights physical groups equally; MAE and bias weight rows equally.

The same unchanged physical formulation must be matched across conditions. Standards used for correction are excluded from held-out validation.

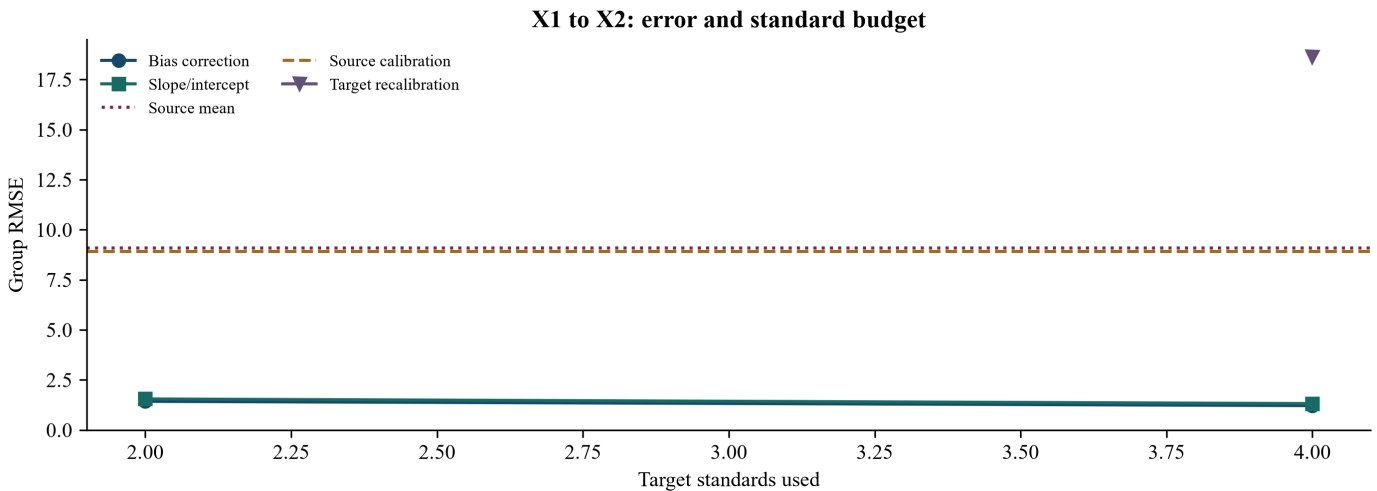

Errors use Acetaminophen content (% w/w), on the recorded common held-out cohort. Lines join evaluated budgets only; they do not imply a continuous optimum. Unadapted source baselines use zero target standards.

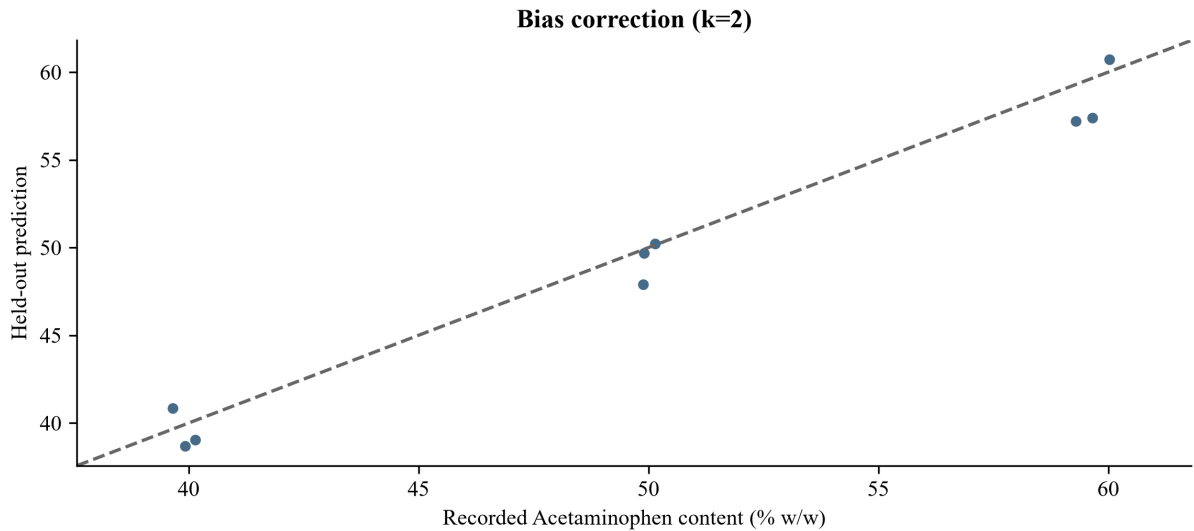

All 9 finite pairs displayed. Dashed line is equality. Acetaminophen content (% w/w). Scope: common comparison cohort; exactly the 9 headline-score rows. All 9 available pairs remain in extended predictions; summary_cohort.csv records membership.